




# MPI Short<sup>®</sup> System

We implant a better future

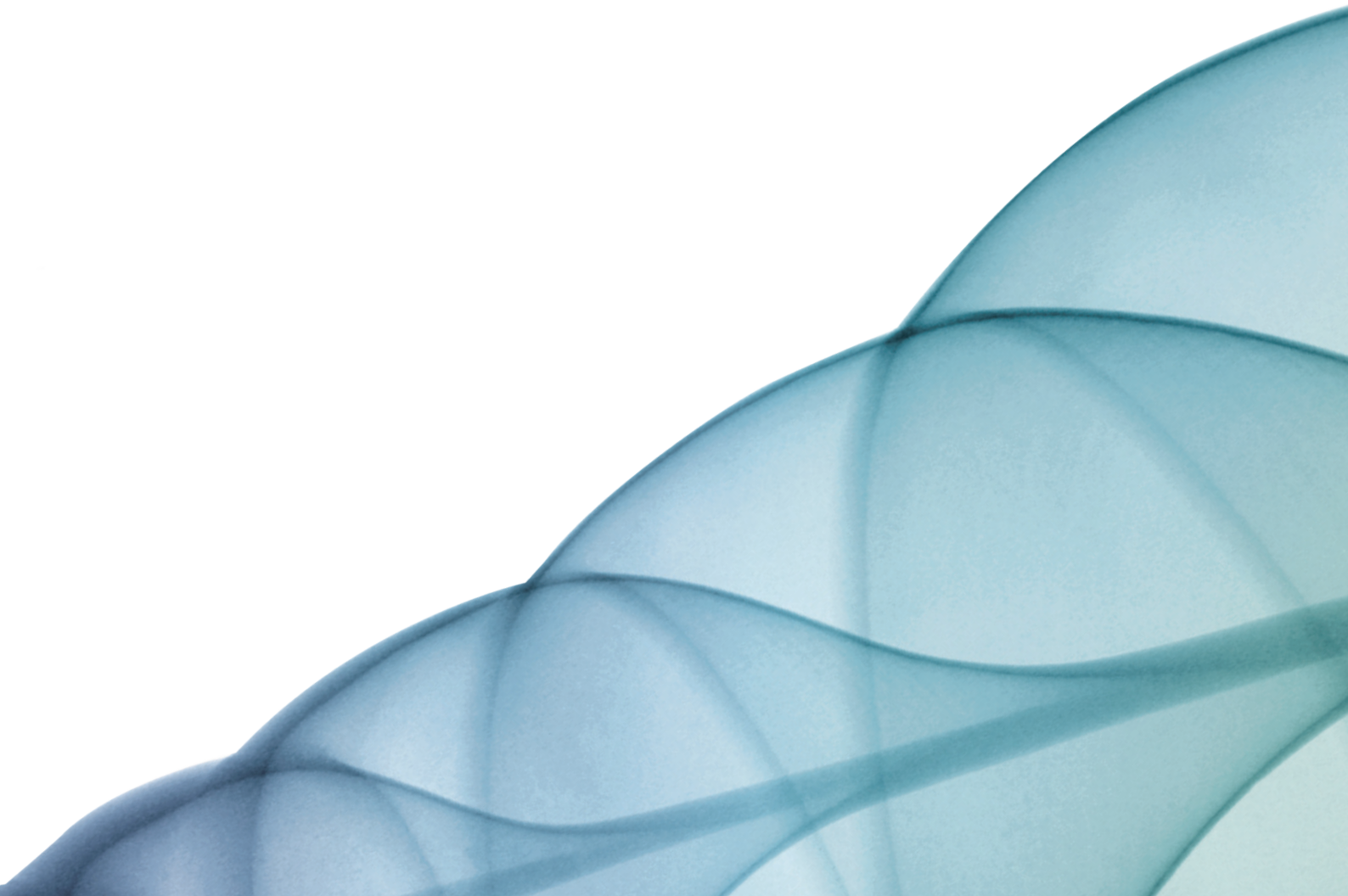

## Advantages MPI Short®

The MPI Short® dental implant has a straight and partially conical body with an active thread profile that provides good primary stability. The length and characteristics of the MPI Short® dental implants are designed to offer an optimal solution for compromised cases, in which the bone height is insufficient to accommodate standard length dental implants.

### Coronal Design

Designed with a micro thread to reduce the reabsorption of the vertical crestal bone. This characteristic improves the load distribution and its fatigue resistance.

The design allows to obtain higher values of initial stability due to the osteotome effect generated during its placement. The platform switching is guaranteed by the machined level of 0,4 mm, which creates a biological width that helps preserve bone levels.

### Thread Design

The double entrance of the threads enables an easy and continuous insertion of the dental implant. The active and self-tapping threads increasing in the apical direction condense the bone gradually, obtaining high initial stability quotients in any type of bone densities.

### Apical Design

Its double helical cutting sides, as well as its active apex allow an axial condensation of the bone, achieving a gradual osteotome effect.

INTERNAL  
CONNECTION  
CONICAL SEAL

### Conical Seal Design

The double internal hexagon permits a fast and adjusted insertion of the prosthetic components, achieving a completely hermetic mechanical closing and tightness.

The mono-block behaviour prevents bacterial filtrations and strengthens the design avoiding the physical fatigue when they work as a set.

### Leeway to Manoeuvre and React

MPI has designed a machined surface in the coronal area of the implant in order to be able to react to bacterial colonizations in initial stages.

Decontamination processes are more effective and faster and a minimization of the bone loss is granted.

### Surface Treatment Integratec®


Double surface treatment. Its optimal macro- and microstructure favours the adherence of osteoblasts and increases the bone-implant contact.

### High Performance Titanium


The Grade IV Titanium used for the MPI Short® dental implants is cold formed. This treatment of the material improves its hardness by 20%, thus increasing the resistance of the implant to higher insertion torques. Material of German origin.

## Internal Connection Conical Seal

### Product features

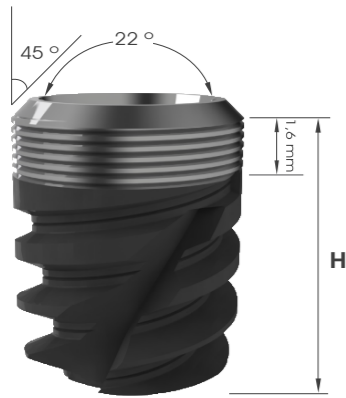

- The microthreads in the coronal area help to reduce the reabsorption of the crestal bone.
- The height of the active threads allow an optimal anchoring in reduced bone situations.
- Two helical cutting sides condense the bone gradually, obtaining high initial stability quotients in all types of bone densities.
- Maximum resistance of the implant and its connection due to the use of cold-formed Titanium Grade IV.
- The placement is made with the hand piece connector which goes directly into the implant.
- Suitable for restorations in the posterior area with reduced bone heights.

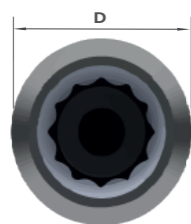

		D - Diameter (mm)
H - Length (mm)	ISCM	Ø 5
	6	ISCM56
	7	ISCM57
	8,5	ISCM58



### Cover Screw

- Made of Titanium Grade V with an anodizing treatment.
- Use with screwdriver hexagonal manual Ø 1,2 mm (SDHM 121 / SDHM 122).

Ø 5
SCM595

## Drilling sequence ISCM

		Diameter of implant (mm)
		Ø 5
	Lance drill	
	Drill Ø 2,0 mm	
	Drill Ø 3,0 mm	
	Drill Ø 4,0 mm	
	Drill Ø 4,8 mm	
	Drill Ø 4,8 mm straight	

DRILL IDENTIFICATION AND BATCH CODE

10 mm

8,5 mm

7 mm

6 mm

0 mm

0,3 mm

BONE TYPE IV	BONE TYPE III	BONE TYPE II	BONE TYPE I
--------------	---------------	--------------	-------------

Prosthetic Components

**Internal**  
Connection  
Conical Seal



## Healing Abutment Straight

- Made of Titanium Grade V.
- Maximum torque recommended 10 N-cm.
- Use with screwdriver 1,2 mm (SDHM121 / SDHM122).

	Ø 5
2 mm	SCM505
3 mm	SCM506
4 mm	SCM507



## Cr-Co Base Castable Abutment

- Maximum torque recommended 30 N-cm.
- For direct implant restorations.
- For cemented prosthesis.
- All non-rotating castable abutments are in black color and the rotating castable abutments in white color.
- Served with abutment screw (SCM536).

	HEXAGONAL	CYLINDRICAL
Ø 5	SCM576	SCM578



## Impression Coping Open Tray

- Made of stainless steel.
- Maximum torque recommended 10 N-cm.
- Use with screwdriver 1,2 mm (SDHM121 / SDHM122).



## Ti-Base

- Manufactured in Titanium Grade V and Titanium Nitride (TiN) coating.
- Served with the abutment screw (SCM536).

	HEXAGONAL	CYLINDRICAL
Ø 5	SCM601	SCM603



## Implant Analog

- Made of stainless steel.

Ø 5
EXC750



## Scanbody

- Made of Titanium Grade V.
- Surface treatment to avoids reflection.
- Abutment screw (SCM536) included.

Ø 5
SCM611

## Abutment Screw

- Made of Titanium Grade V.




Ø 5
SCM536

# General Information

## Packaging Opening Instructions

Designed for easy identification of each implant. According to the color of the cardboard box, we can find:

- Blue packaging: Implant MPI Short® CM

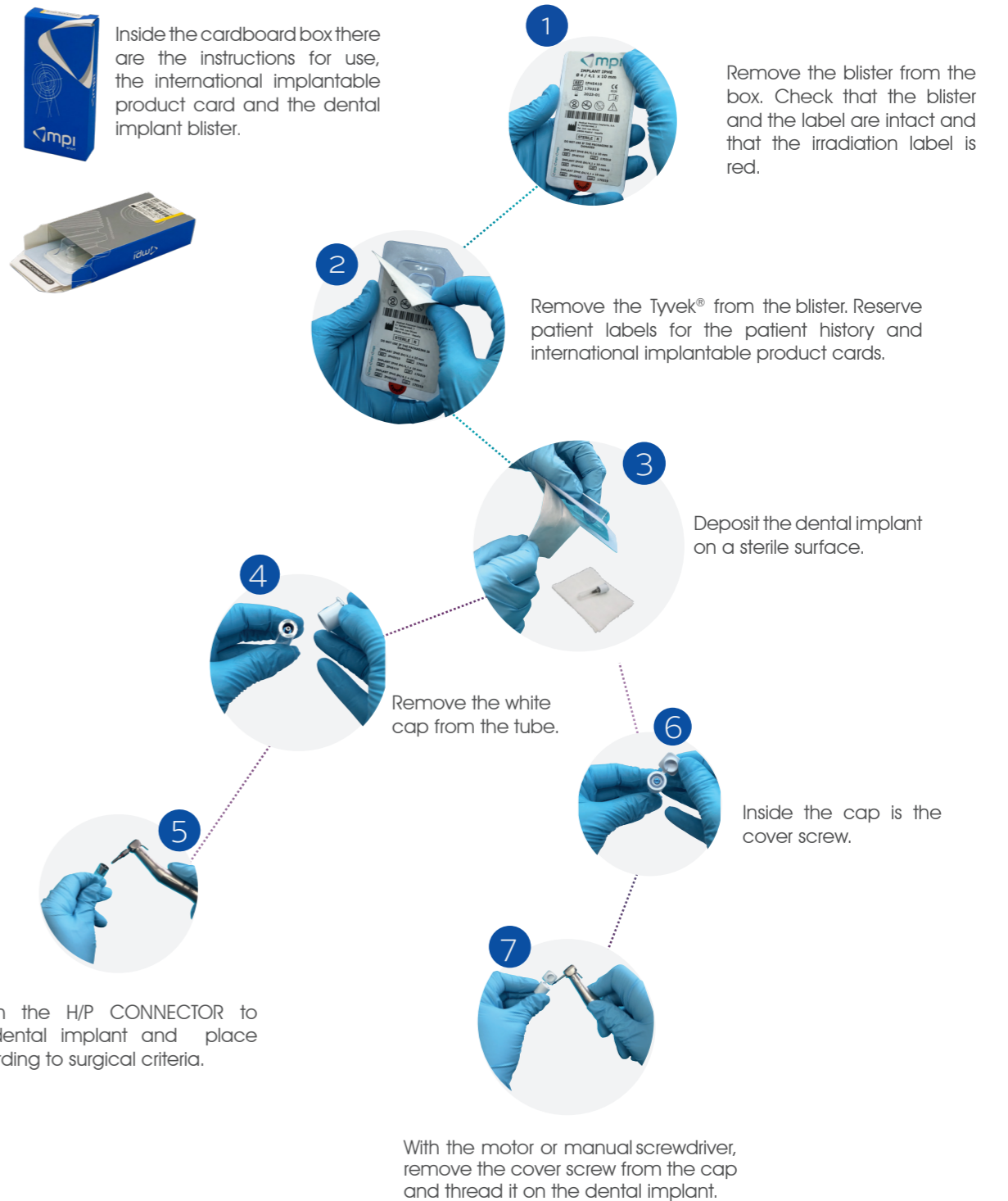

# International Implantable Product Card

# Label Symbology

MPI will deliver with each of its implantable products, an international implantable product card.

This card contains the following pre-printed information:

- MD: Medical device to which the card belongs
- UDI: Unique product identifier.

The user must stick a patient label (included on the product label) on each of the two boxes identified with:

- LOT / UDI-DI.

Front view of the International Implantable Product Card. It features four main sections: MD (Medical Device), LOT (Lot), UDI-DI (UDI-Device Identifier), and UDI (Unique Device Identifier). Each section has a corresponding label and a blank space for information. At the bottom, there is a logo for Medical Precision Implants, S.A. and the company address: C/EMPLEO, 21 P. I. LOS OLIVOS 28906 GETAFE, MADRID, SPAIN. Multilingual text at the bottom reads: 'es Tarjeta de Implante para Paciente/ en Patient Implant Card/ fr Carte d'implantation pour le Patient/ it Scheda di Implanto per Paziente/ pt Cartão de Implante para Paciente'.

Back view of the International Implantable Product Card. It features four main sections: MD (Medical Device), LOT (Lot), UDI-DI (UDI-Device Identifier), and UDI (Unique Device Identifier). Each section has a corresponding label and a blank space for information. At the bottom, there is a logo for Medical Precision Implants, S.A. and the company address: C/EMPLEO, 21 P. I. LOS OLIVOS 28906 GETAFE, MADRID, SPAIN. Multilingual text at the bottom reads: 'es Tarjeta de Implante para MPI/ en MPI Implant Card/ fr Carte d'implantation pour MPI/ it Scheda di Implanto per MPI/ pt Cartão de Implante para MPI'.

On the back of the card, the user must fill in the following information with the patient's information:

- Patient name (Patient ID)
- Record number (Patient number)
- Implementation date
- Name of the medical care center.


Patient label for the International Implantable Product Card. It features a header with 'Medical Precision Implants, S.A. International Implant Card'. Below the header are four lines for patient information, each with an icon: a person with a question mark for name, a number '31' for record number, a person with a plus sign for implementation date, and a person with a plus sign for medical care center. At the bottom, there is a logo for Medical Precision Implants, S.A. and the website address: www.mpimplants.com. Multilingual text at the bottom reads: 'es para Paciente/ en for the Patient/ fr pour le Patient/ it per le Paziente/ pt para Paciente'.

Patient label for the International Implantable Product Card. It features a header with 'Medical Precision Implants, S.A. International Implant Card'. Below the header are four lines for patient information, each with an icon: a person with a question mark for name, a number '31' for record number, a person with a plus sign for implementation date, and a person with a plus sign for medical care center. At the bottom, there is a logo for Medical Precision Implants, S.A. and the website address: www.mpimplants.com. Multilingual text at the bottom reads: 'es para MPI/ en for MPI/ fr pour MPI/ it per MPI/ pt para MPI'.

The international implantable product card marked with:

- "for the Patient", must be given to the patient and
- "for MPI", must be returned to MPI by the medical care center.

Label:

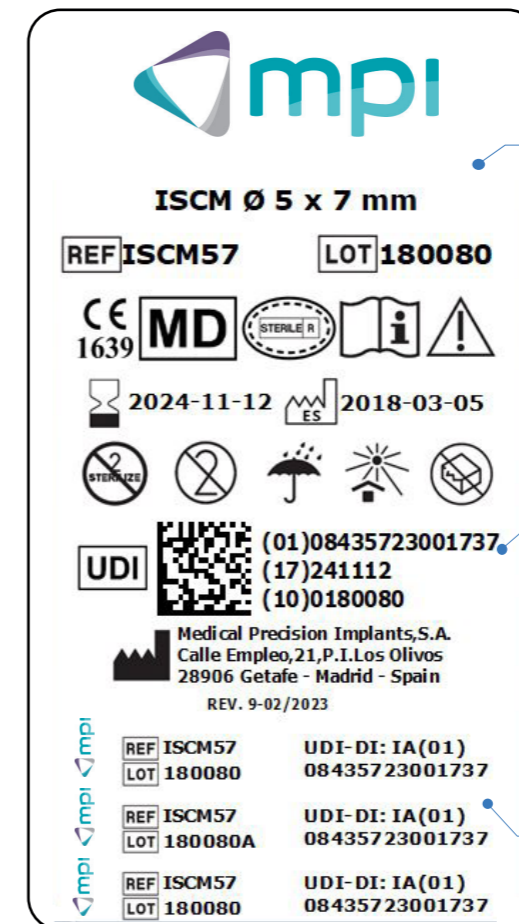

Dental implant description and measurements are indicated:


- ISCM: MPI Short® System

UDI - DataMatrix  
(01) UDI-DI  
(17) Expiration date (YYMMDD)  
(10) Batch number

Patient labels indicating reference, batch code and UDI-DI.

Symbology:

	Manufacturer.		Caution.
	Country of manufacture with the two-letter country code defined in the Standard ISO 3166-1. The date of manufacture can be added adjacent to the symbol.		UDI-DI Data matrix: (01) UDI-DI / (17) Expiration date (YYMMDD) / (10) Batch number.
	CE symbol with the notified body number 1639.		Use by date.
	Batch code.		Sterilized using irradiation. Single sterile barrier system with protective packaging inside.
	Catalogue number.		Unique Device Identification.
	Keep away from sunlight.		Medical Device.
	Do not re-use.		Text language.
	Do not re-sterilize.		Do not use if package is damaged and consult instructions for use.
	Keep dry.		
	Consult instructions for use.		

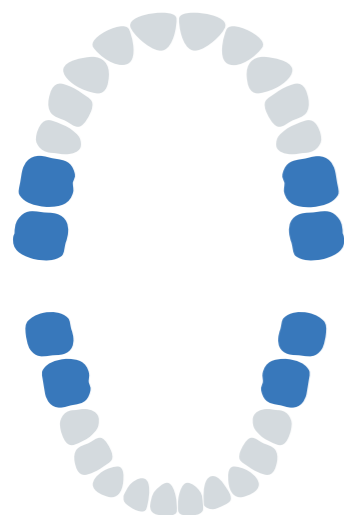

**Implants diameters:** Wide platform.

- MPI Short® CM Ø 5

**Instructions for use:**

For fixed, multiple restorations.

**Recommended position:**

Molars in both jaws.

**Recommendations:**

- Recommended for multiple restorations: the union of several implants allows an increase in resistance to masticatory forces and a distribution of loads on the bone surface.
- Maximum insertion torque: 80 N·cm.
- Insertion rotation: 30 rpm.
- Immediate loading is not recommended, always place in two-phase surgery.
- Under-drilling recommended to obtain an osteotome effect.

**Quality**

To ensure the most demanding quality standards, we have the help of our metrology laboratory equipped with qualified personnel and state-of-the-art measuring instruments. Thanks to demanding quality controls, we have achieved a success rate of over 99%. The “Zero defects” policy is our fundamental basis for the quality of our products and the satisfaction of our customers.

MPI is certified with:

- ISO13485:2016
- UNE166002:2021

**Orders**

Telephone: +34 91 684 60 63

E-mail: pedidos@mpimplants.com

**Shipments and verification of the products**

For orders over 150€, the products will be sent free of charge.

Once the order is received, the user should check the order and products.

If there should be any deviation in the order, please contact us as soon as possible.

**Return policy**

The maximum return period is 15 days from the date of issuance of the delivery note. Changes will not be accepted once the material has been opened and used or its packaging is not in perfect condition.

Returns of unsuccessful materials will not be accepted after 2 months after the cause.

**Payment Methods**

Direct debit / cash on delivery (additional cost).

**Warranty**

We guarantee the replacement of our products as long as they are combined with original MPI components.





For more information:  
☎ +34 916 846 063  
[www.mpimplants.com](http://www.mpimplants.com)